

Data sheet


Cooled Incubator (IL) ILW 115 Smart



The photo above is for reference only, may show additional options not included in standard equipment. The real appearance, particularly color and structure of the material may differ from the ones presented in the photo.

Advantages of the SMART controller:

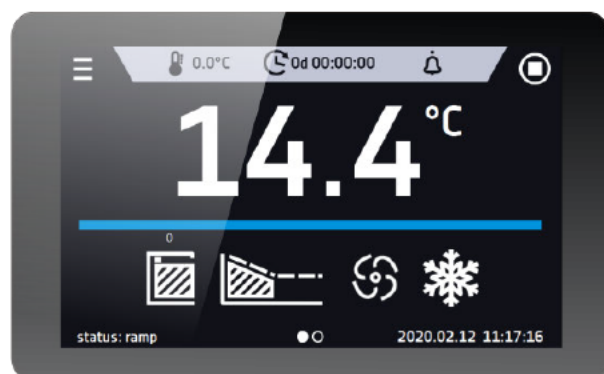- 4,3", clear, full colour touch screen
- LAN, USB ports for data transfer
- multi-segment time and temperature programs
- visual and sound alarm
- internal memory for programs and data storage
- event registry
- user manual for direct download
- Quick change of program parameters
- Alarm Bar
- operating with gloves on



Smart - preview screen

TECHNICAL DATA

air convection	forced
chamber capacity [l]	112
working capacity [l]	112
controller	microprocessor PID
display	4,3 full colour touch screen

TEMPERATURE


temperature range [°C]	0...+70 / -10...+70 (option)
temperature resolution every ... [°C]	0,1
temperature fluctuation at 37°C [±/°C]*	0,2
temperature variation at 37°C [±/°C]*	0,3
over temperature protection	class 2.0 to DIN 12880 / class 3.3 (option)

CHAMBER

door type	double(5) / door with viewing window (option)
interior	
Smart	DIN 1.4301
IG Smart	DIN 1.4301
housing	
Smart	powder coated sheet
IG Smart	stainless steel linen finish

overall dims [mm] /1/

width A	660
height B	1140
depth C	720
internal dims [mm]	
width D	460
height E	540
depth F	450

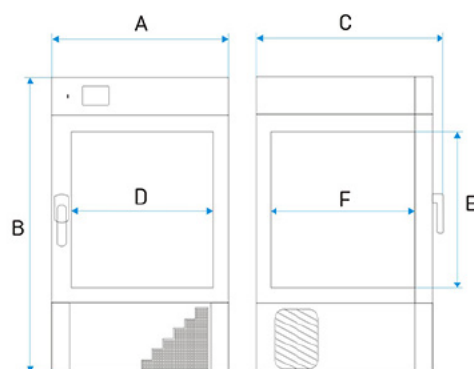

shelves (standard max)	2 7
max shelf workload [kg] /2/	25
- reinforced shelf version (PW) [kg] /3/	50
max unit workload [kg]	60
- reinforced unit version (W) [kg] /4/	120
weight [kg]	90

ELECTRICAL PARAMETERS

voltage**	230V 50-60Hz
nominal power [W]	500
refrigerant	1234ze / GWP=1
warranty	24 months
manufacturer	POL-EKO-APARATURA
txt_opis pod tabelą	<p>all the above technical data refer to standard units (without optional accessories)</p> <p>* - fluctuation measured in centre of the chamber; in space, variation (K) calculated for chamber as: $K = \pm (T \text{ average max.} - T \text{ average min.}) / 2$</p> <p>** - other power supplies on request</p> <p>1 - depth doesn't include 50 mm of power cable 2 - on uniformly loaded surface 3 - reinforced shelf 4 - reinforced version 5 - additional internal glass door</p>

all the above technical data refer to standard units (without optional accessories)

* - fluctuation measured in centre of the chamber; in space, variation (K) calculated for chamber as:

$K = \pm (T \text{ average max.} - T \text{ average min.}) / 2$

** - other power supplies on request

1 - depth doesn't include 50 mm of power cable

2 - on uniformly loaded surface

3 - reinforced shelf

4 - reinforced version

5 - additional internal glass door

OPTIONS AND ACCESSORIES



Order number: */A

door with viewing window



INOX

Order number: */P INOX

Stainless steel wire shelf INOX



Order number: */PP

Perforated shelf



Order number: */PW

Reinforced shelf



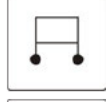
Order number: KUW GN**

Stainless steel cuvettes



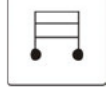
Order number: QLK*

Wheels











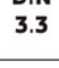



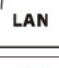
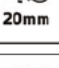


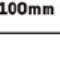
Order number: */S or */S INOX

Table with wheels



Order number: */ST or */S INOX

Base on castors

	Order number: *PLUS	Automatic defrosting function
	Order number: *T	Low temperature version
	Order number: */W	Reinforced version
	Order number: OWW/OWW LED	Interior lighting
	Order number: */FOT	Photoperiodic system
	Order number: GNZ	Internal socket
	Order number LabDesk	LabDesk software
	Order number: BRT*/L or IQ/OQ/PQ	Calibration and IQ, OQ, PQ qualification
	Order number: */3.3	Over temperature protection 3.3 class according to DIN 12880
	Order number: BPP 12	Battery backup for display
	Order number: PORT ALARM	Dry alarm contact
	Order number: USBK	USB cable
	Order number: LANK	LAN cable
	Order number: OCZ/20	Non-standard access port 20 mm
	Order number: OCZ/30	Additional access port 30 mm
	Order number: OCZ/60	Non-standard access port 60 mm
	Order number: OCZ/100	Non-standard access port 100 mm

ALEPH

Part Science · Part Art



Order number: KD

Access control